



## Anonymizer – Batch processing

It anonymizes the set {CTset, StructureSET, Plan, Dose} preserving the logic relationship between objects.


The anonymization is done using the MATLAB function *dicomanon* (rev R2010a) which anonymizes as per stipulated on [PS 3.15-2008 Table E.1-1](#). We are not able to test full conformity of this MATLAB tool to DICOM de-identification standards, thus we do not take any responsibility for the use of this tool.

For keeping track of some possible useful details, the Patient Age is preserved declaring a fake Birth Date. Patient sex is also preserved.

How to use:

Create a “root” folder on which each set {CTset, StructureSET, Plan, Dose} will be stored on a different subfolder of that root folder.

Export each set {CTset, StructureSET, Plan, Dose} to a single sub-folder. Only 1 plan, 1 structure set and 1 Dose are allowed to be on the folder with the CT image set.

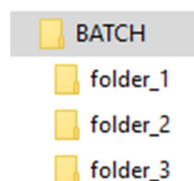

The software does not check for consistency of this set, thus is your responsibility to export a consistent set, i.e. the structure set defined on the CT set, the plan defined on the structure set and the dose corresponding to this plan.


Set the new generic PatientID, PatientFirstName, PatientLastName and the name of the Hospital of origin of the set {CTset, StructureSET, Plan, Dose}.

Set the Initial Index #: The PatientID and PatientLastName will have the Index # added at the end, in sequential order starting from Initial Index #: NEWID005, NEWID006, etc.

Once edition completed and checked, press the “Check data” to enable the “Anonymize” button that has to be pressed to launch the process.

The dialog box asks for selection of the root folder.

The result is a series of new folders **.\ANON\_###** (created on the root folder) containing the new anonymized DICOM objects.

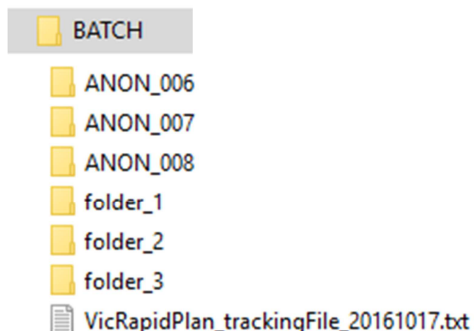



A text file named “VicRapidPlan\_trackingFile\_yyyymmdd.txt” is created on the root folder.) containing the following information: Original PatientID, New PatientID, Original PlanID, Hospital of origin. The first line have the headers of the data and following lines the data itself. Its purpose is to facilitate managing the plan flow into models. Be aware it contains the Original Patient ID.

The original DICOM files are not changed.

The new set makes reference to a machine named “VicRapidPlan” which is a copy of the Varian testing machine “EclipseCAP\_TB”, a TrueBeam machine that contains many different beam types. Very likely most of the centers have already this CAP machine declared into their ARIA systems, if not, the XML descriptor of this machine is provided (VicRapidPlan\_TB.xml). You can import this machine for easing the managing of plans with energies not matching the ones declared on your system.